




Mirage Micro™

NASAL MASK

PROVEN SEAL, PERFECT FIT

The latest generation nasal mask from ResMed provides personalised fit, seal and comfort for more users. Working with the proven Mirage™ dual-wall cushion,

the mask's unique MicroFit dial fine tunes the fit to suit each individual user. Intuitive to use by sleep professionals and users, Mirage Micro is the new benchmark in nasal masks.

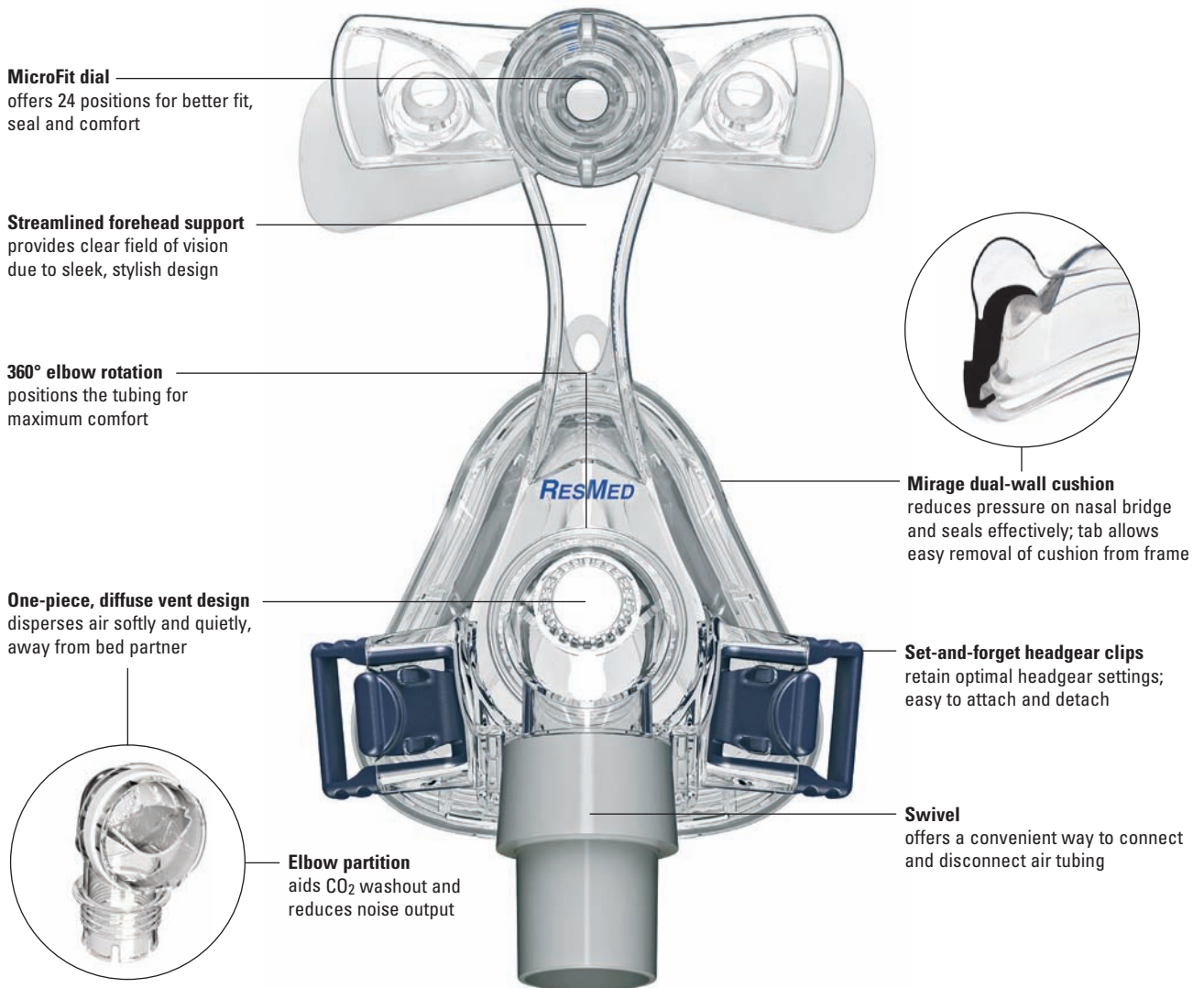

MicroFit dial
offers 24 positions for better fit, seal and comfort

Streamlined forehead support
provides clear field of vision due to sleek, stylish design

360° elbow rotation
positions the tubing for maximum comfort

One-piece, diffuse vent design
disperses air softly and quietly, away from bed partner



Elbow partition
aids CO₂ washout and reduces noise output



Mirage dual-wall cushion
reduces pressure on nasal bridge and seals effectively; tab allows easy removal of cushion from frame

Set-and-forget headgear clips
retain optimal headgear settings; easy to attach and detach

Swivel
offers a convenient way to connect and disconnect air tubing



TECHNOLOGY HIGHLIGHT

MicroFit dial

A choice of 24 positions at your fingertips—providing a comfortable, effective fit for more users

- **Adjustable:** can easily be adjusted while on the user's face
- **Precise:** each turn of the dial adjusts the forehead support in small increments
- **Progressive:** moves in towards the face, allowing the user to sense when the fit is most comfortable
- **Works in combination with the cushion:** as the dial adjusts the angle of the mask to the face, the dual-wall cushion adapts to enhance the seal without the need for overtightening the headgear

USER BENEFITS

Intuitively easy

- MicroFit dial is easily adjusted with one hand even while lying down
- Vent cover can be removed for easy and hygienic cleaning
- Intuitive mask design for easy handling, cleaning, assembly and disassembly

Personalised fit and comfort

- MicroFit dial fine tunes the fit to suit each user's profile
- Mirage dual-wall cushion is available in five sizes, including new Small size, to fit more users
- Deep mask frame design accommodates a greater range of nasal structures

Quiet, sleek and stylish

- Array of small, built-in vents diffuse air quietly
- Streamlined design feels light on the face

PROFESSIONAL BENEFITS

First-time fit for more users

- Combination of MicroFit dial and dual-wall cushion achieves a perfect fit for a range of facial structures—first time, every time
- Five cushion sizes including new Small size (cushion and frame) fit a wider range of users

Customer satisfaction and compliance

- Increased therapy acceptance through first-time fit
- Personalised fit, comfort and seal encourage long-term compliance
- Mask parts designed for fast and intuitive assembly and disassembly
- Comprehensive user support package including travel bag and user guide with mask fitting and care instructions

Cost and time efficiencies

- Reduced fitting and training time: dial and cushion combination seals quickly
- Reduced callbacks: ease of use and fitting
- Two frame sizes fit five cushion sizes

PRODUCT CODES

Mirage Micro Nasal Mask

Europe¹

Small	16338
Medium	16339
Large	16340
Large Wide	16341
X-Large	16342

Europe²

Small	16343
Medium	16344
Large	16345
Large Wide	16346
X-Large	16347

Europe³

Small	16348
Medium	16349
Large	16350
Large Wide	16351
X-Large	16352

Asia Pacific

Small	16368
Medium	16369
Large	16370
Large Wide	16371
X-Large	16372

¹ English, French, German, Spanish, Italian, Portuguese and Dutch

² English, Swedish, Norwegian, Danish, Finnish, Hungarian, Polish and Greek

³ English, Russian, Turkish, Czech