



Ultra Mirage™ II Nasal Mask


The all-round performer

ResMed's Ultra Mirage™ II Nasal Mask is suitable for the majority of patients on positive airway pressure—using ResMed's CPAP, bilevel, or AutoSet™ devices.

Available in four cushion sizes (Standard, Large, Shallow, and Shallow-wide), the Ultra Mirage II will fit most patients—first time, every time.


The Ultra Mirage II continues to set the benchmark for nasal masks in terms of seal, comfort, and quietness—valuable for effective therapy.

It incorporates ResMed's renowned Mirage dual-wall cushion technology and a unique air vent design for ultra-quiet performance.

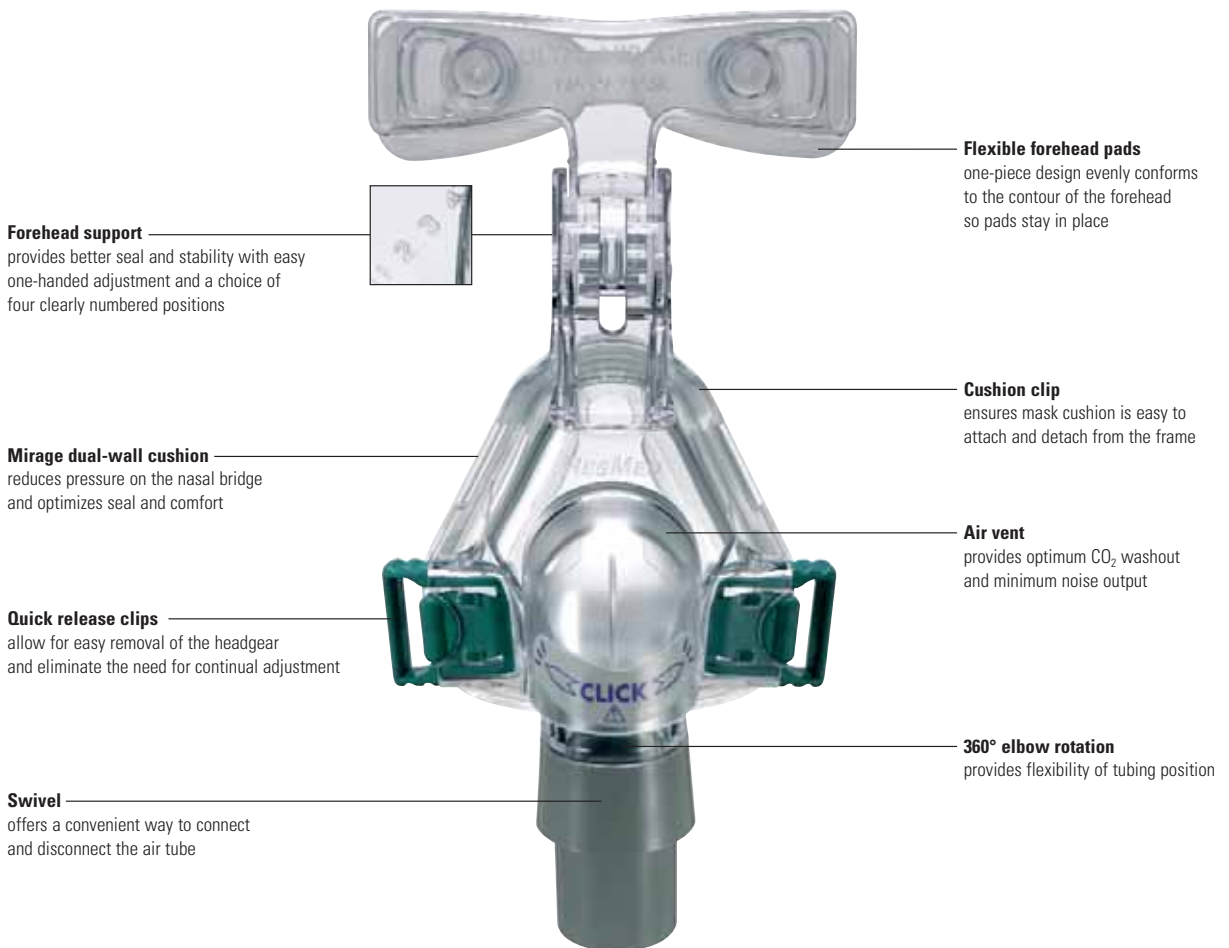

TECHNOLOGY HIGHLIGHT— UNIQUE AIR VENT

- Unique air vent design minimizes airflow intensity and noise.
- Air from the vent is directed down the tube, away from both patient and partner.
- Air vent specially designed to provide adequate flow to wash out exhaled CO₂.

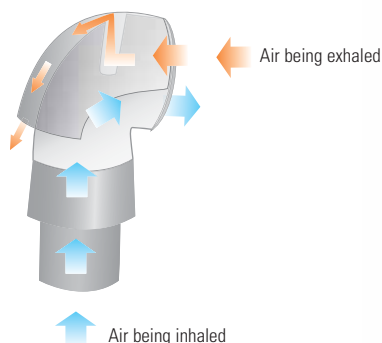

Adjustable forehead support

- Provides additional support especially for those patients on higher pressures and/or bilevel therapy.
- Offers a choice of four numbered positions for better seal and stability.

Swivel

- Provides a convenient way to connect and disconnect the air tube, without removing the mask.
- Allows flexibility for patient movement.

Rotating elbow

- 360° elbow rotation allows patients to choose air tube position for maximum comfort.



Adjustable forehead support—four clearly numbered positions

Secure and comfortable headgear

- Headgear clips allow quick attachment/detachment of headgear—no further adjustments required after initial fitting.
- Four points of attachment to mask frame for enhanced stability.

Quick fitting

- Mask is packaged fully assembled with headgear.
- Adjustment of forehead support and headgear is all that is required on initial setup.



Rotating elbow—flexibility of air tube position

ULTRA MIRAGE II NASAL MASK SYSTEM PRODUCT CODES

North & Latin America

Standard	16548
Large	16549
Shallow	16550
Shallow-wide	16577

UK

Standard	16321
Large	16322
Shallow	16323
Shallow-wide	16324

Europe (except UK) and Asia

Standard	16545
Large	16546
Shallow	16547
Shallow-wide	16578

Australia & New Zealand

Standard	16579
Large	16580
Shallow	16581
Shallow-wide	16582