

**ResMed**

**AirMini™**

## The revolution you've been waiting for

Introducing ResMed AirMini, the world's smallest CPAP device ideal for patients looking for compact and discreet therapy.

Featuring a unique venting system and waterless humidification, it combines our latest technologies in a sleek, pocket-sized design that can be used at home or while patients are away.

To meet patients' lifestyle needs, satisfy comfort expectations and boost business benefits, start small with AirMini.



## Small. Smart. Proven.

AirMini leverages the same proven algorithms as our market-leading bedside devices to help treat patients weighing 30 kg or more with obstructive sleep apnea (OSA). It also includes ResMed AutoSet™ for Her, a mode which has been specifically designed to treat women with mild to moderate OSA.

With Bluetooth® connectivity and a dedicated app, AirMini helps patients and clinicians set up for success and stay engaged with ongoing therapy.

### Key features


- ResMed HumidX™ waterless humidification minimises heat and moisture loss for comfort throughout the night.
- ResMed ActiveAir™ advanced venting enables our latest therapy technologies to be delivered in a compact system.
- Bluetooth connectivity to smart devices and the dedicated AirMini App offer user-friendly setup and use.



AirMini Travel Bag  
**38840**



AirMini Travel Case  
**38841**

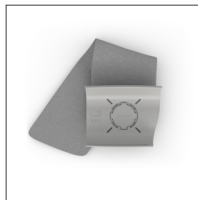


AirMini Mount System  
**38842**



HumidX **38808** (50 pk)  
HumidX **38809** (3 pk)  
HumidX **38810** (6 pk)



HumidX Plus **38811** (50 pk)  
HumidX Plus **38812** (3 pk)  
HumidX Plus **38813** (6 pk)

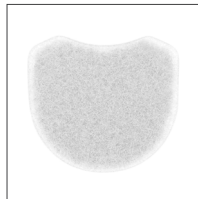

AirMini Filter  
**38837** (50 pk)  
AirMini Filter Hypo  
**38859** (50 pk)

### Product codes

AirMini Autoset APAC	<b>38140</b>
N20 Setup Pack	<b>38823</b>
P10 for AirMini	<b>38824</b>
F20 Setup Pack	<b>38825</b>

### Technical specifications

#### AirMini device

Dimensions:	133mm (W) x 84mm (D) x 52mm (H)
Weight:	300g
Housing construction:	Flame retardant engineering thermoplastic

#### Modes

Operating pressure range	4 to 20 cm H <sub>2</sub> O (4 to 20 hPa) in AutoSet, AutoSet for Her and CPAP modes
--------------------------	--

#### AirMini power

Standby power consumption:	1W
----------------------------	----

#### Sound

Pressure level measured according to ISO 80601-2-70:2015 (CPAP mode) is 29 dBA with uncertainty of 2 dBA.
---

#### Power supply unit

AC input range:	100–240V, 50–60Hz 115V, 400Hz (nominal for aircraft use)
DC output:	24V, 0.83A
Typical power consumption:	6.3W
Peak power consumption:	27W

#### Environmental conditions

Operating temperature:	+5°C to +35°C
Operating humidity:	10 to 95% relative humidity, non-condensing
Operating altitude:	Sea level to 2,591 m; air pressure range 1013 hPa to 738 hPa
Storage and transport temperature:	-25°C to + 70°C
Storage and transport humidity:	5 to 95% relative humidity, non-condensing

#### Aircraft use

ResMed confirms that this machine meets the Federal Aviation Administration (FAA) requirements (RTCA/DO-160, section 21, category M) for all phases of air travel.
--

#### Electromagnetic compatibility

The AirMini complies with all applicable electromagnetic compatibility requirements (EMC) according to IEC60601-1-2:2014, for residential, commercial and light industry environments.
--

#### IEC 60601-1:2014 (Edition 3.1) classification

Class II (double insulation), Type BF, Ingress protection IP22.
---